verilux®

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Portable LED Light Therapy Lamp

Model: VT54

Welcome to the Verilux Family!

We're so excited that you have purchased this HappyLight® Therapy Lamp. Why? Because light therapy has brightened the lives of millions of people, and you now hold a product in your hands that has the power to brighten your life, too. Pretty amazing stuff!

We want to make sure that you get the most out of your HappyLight® Therapy Lamp. That's why we stand behind our products and have a three-year limited warranty on this lamp. It's also why we have a super supportive Customer Happiness Team and resources on our website to help you create a light therapy routine that's completely personalized.

Call us at 1-800-786-6850 or visit us at verilux.com if you have any questions as you get started. We're here to help, and we love hearing from customers like you.

Have a bright day!

The Verilux Customer Service Team

Contents

FCC Notice.....	4
Important Safeguards.....	5
Setup.....	6
Device Pretest.....	6
How to Use.....	7-11
Controls.....	7
Stand and Positioning.....	8-9
Customization.....	10
Establishing a Daily Routine.....	11
Features and Benefits.....	11-12
Care and Cleaning.....	12
Technical Specifications.....	12
Troubleshooting.....	13
Warranty.....	13-14

FCC Notice

Model: HappyLight® Mini Portable LED Light Therapy Lamp (VT54)

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Verilux
2803 S. Yale St.
Santa Ana, CA 92704
verilux.com

The product is protected against unauthorized use. DO NOT modify it without the authorization of the manufacturer, or the user's authority to operate could be voided.

Note: This equipment has been tested and found to comply with the limits for a Class B[†] digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

[†] This Class B digital apparatus complies with Canadian ICES-005.

READ ALL INSTRUCTIONS BEFORE USING

Important Safeguards

⚠ DANGER:

- To avoid electrocution, do not operate this device near water.

⚠ WARNING:

- Close supervision is recommended when this lamp is used by or near children, or others who may not understand the need for these precautions.
- To prevent risk of shock or personal injury when cleaning the lamp, make sure you have turned off and unplugged the lamp and allowed sufficient time for the bulb to cool.
- Do not shorten or cut the power cord.
- Certain medications (e.g., anti-depressants) may induce light sensitivity which may result in discomfort or a level of hyperactivity. Consult with your health care provider before using this product.
- If you have a pre-existing eye condition or are at risk of developing Age-Related Macular Degeneration or other eye conditions, consult with your healthcare provider before using this product.
- Some people experience light sensitivity to bright light or upon over-exposure to bright light. Discontinue use if you experience any discomfort.

⚠ CAUTION:

- This lamp is for indoor use on a hard surface (e.g., a table or desk) only.
- Never cover the lamp or place anything on top of it when it is in operation.
- Do not operate this lamp if it has been damaged in any way.
- Do not disassemble. There are no user-serviceable parts in this lamp.
- For best compatibility, it's recommended to use the power cord and adapter supplied with the lamp. Unit may not be compatible with all other power adapters, cords, or connectors.
- This lamp is not for use with light dimmers, timers, motion detectors, or voltage transformers.
- Use a quality surge suppressor on all electronic equipment. Voltage variations and spikes can damage electronic components in any system.
A quality suppressor can eliminate the vast majority of failures attributed to surges. Suppressors may be purchased at electronics stores.

SAVE THESE INSTRUCTIONS

Setup

Plug In

- Place lamp on a flat, level surface.
- Insert the USB-C connector of the cable into the USB-C port - as shown on the back of the lamp.
- Insert the USB-A connector into the adapter and plug into a standard 120V AC electrical outlet.
- Turn on the lamp using the power button.

Note: HappyLight® Mini must be plugged into the included power adapter or another 3A power source to achieve 10,000 LUX. Standard wall adapters with lower power ratings (<2A) will not provide the full 10,000 LUX illumination necessary to achieve the benefits of light therapy.

Device Pretest

To help determine your degree of sensitivity to light, we recommend that you perform this pretest upon first use of your HappyLight® Mini Therapy Lamp.

Step 1 - Turn on the lamp and position it about 18 inches from your face (remember, don't look directly into the lamp).


Step 2 - Operate the lamp for approximately 10 minutes and turn it off. If the light has caused any discomfort, stop using the lamp and consult a doctor. If you have no discomfort, you may consider gradually increasing the brightness of the light, decreasing the distance from the lamp, and/or increasing your usage time. When you first start using the lamp, you should monitor the degree of any discomfort and accordingly adjust the brightness, distance, and time. Keeping a journal or log of how you feel when using HappyLight® Mini during the first few weeks will help you find the settings, distance, and time that work best for you!

⚠ WARNING: Certain medications (e.g. antidepressants) and diagnosed mental health conditions may induce light sensitivity, which may result in discomfort, hyperactivity, or hypomania. Consult with your healthcare provider before using this product.

⚠ WARNING: If you have a pre-existing eye condition or are at risk of developing Age-Related Macular Degeneration or other eye conditions, consult with your health provider before using this product.


How to Use - Controls

Your HappyLight® comes with touch controls, making your light therapy experience easy and fun!

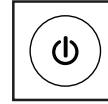
On/Off: The main power button is located on the front of Mini. The light must be turned on before you can change specific settings. When you turn off your light, your selections are automatically saved for the next time you use your light. How cool is that?


There are two controls that allow you to change your settings and customize your HappyLight® experience:

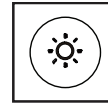
Brightness: Choose between 3 different brightness levels to find the one most comfortable for you*. If the brightest setting is uncomfortable, you can start with a lower setting and work your way up to the full 10,000 lux over multiple sessions. Increasing the brightness when in use will allow you to get your daily dose of healthy, bright light more quickly (which means shorter light therapy sessions!). However, if a lower setting is more comfortable for you, that's OK! Just increase your time with the light and you will still get the same benefits!

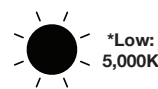
*Low:
5,000K

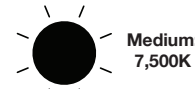Medium:
7,500K

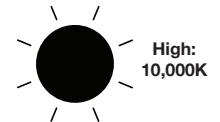High:
10,000K

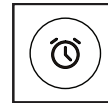Countdown Timer: The countdown timer helps you keep track of the length of your light therapy experience (which makes it pretty convenient). You can set the timer while you do other things (like reading a book, checking email or just enjoying a cup of coffee or tea), and your light will turn off on its own once your light therapy session is complete.

Setting the Timer: Touch the Timer Button to choose between four different timer settings. The lamp will cycle through 15, 30, 45, 60, and off options until you select your preferred time.

15

30

45

60

How to Use - Adjusting Stand

The HappyLight® Mini stand is adjustable up to 90° to ensure you find the perfect angle. Tilt the stand between 15° for a taller angle or adjust the full 90° for a flatter angle.

How to Use - Positioning

Sit close to your HappyLight® Mini with your eyes open. While staying oriented toward the light, feel free to engage in other activities that you would normally do. These can include reading, working at your computer or desk, using your phone or tablet, applying makeup, or even hobbies and crafts! Really, there is no limit to what you can do – get your dose of bright light for the day and do your normal routine.

You can place HappyLight® Mini between 6 and 24 inches from your face. The time of use and distance are correlated - the closer you are, the less time you need in your light therapy session and vice versa when it is placed further away (more time when it is further from you). It is most effective when positioned slightly off-center.

Important Note: Do not look directly into the lamp.

How to Use - Customization

Since the light does need to be oriented towards your face (but not staring into it directly), your eyes need to be open to get the full benefit.

Note: Tinted glasses reduce the amount of light reaching your eyes and would result in the need for a longer session.

We are all different and each person's response to light will vary. It's important to find the distance, time, and brightness that are most comfortable for you. How close you place HappyLight® Mini to you affects the length of time using the light. The closer you are, the higher the brightness and the less time you will need to use HappyLight® Mini each day to experience the benefits!

If you do not experience desired results, or are feeling lethargic, here are a few things you can try:

- Try using a higher brightness level
- Sit closer to the lamp
- Increase the amount of time you are using the lamp

If you experience edginess or feel over-simulated, try some of these options:

- Use a lower brightness level
- Sit farther from the lamp
- Decrease the amount of time you are using the lamp
- Take a break from using the lamp

Try keeping a journal of how you feel when using HappyLight® Mini to find the settings that work best for you!

Establishing a Daily Routine

20 minutes
to
1 hour/day*

Establishing your daily routine with HappyLight® Mini will help you to maximize the benefits of light therapy. Clinical research on bright light therapy consistently demonstrates that 2,500 – 10,000 lux (a measurement of the

intensity of "brightness of light") between 20 minutes and a few hours per day helps with mood and sleep issues and can reverse the symptoms of "winter blues." Consistent use at the same time each day has been shown to have the biggest impact and timing earlier in the day has also been shown to be more helpful for sleep cycles.

Every person has a different sensitivity to light, so it is important to find the settings, distance, and time that are most comfortable for you. After you have completed the Device Pretest, use your HappyLight® regularly around the same time each day. You could choose to do this all in one session, or divide your time into several mini-sessions. Use your HappyLight® Mini daily to improve sleep, boost your mood, increase energy, and enhance your focus.

As you are establishing your routine, if you find yourself feeling sluggish and lethargic, increase the brightness of the lamp, sit closer, and/or spend more time in your session. If you feel a sort of edginess, like you drank too much coffee, it may be time to turn off the lamp.

*Recommended time is a function of brightness and distance. Individual results may vary. Please also refer to the warnings at the bottom of the Device Pretest section on page 6.

Features and Benefits

Your new HappyLight® Mini is designed to give you a light therapy experience that's customizable to fit your unique needs. Here are all the amazing features and benefits you can expect from your new light therapy companion:

- Full spectrum light to brighten up your home, office, or dorm room and cue your body's own energy enhancers to make you more alert, focused, and productive
- 10,000 lux[†] of bright light to improve sleep, boost your mood, increase energy, and enhance your focus

Features and Benefits — continued

- Adjustable brightness with 3 levels that allow you to choose the intensity that you are most comfortable with
- Thanks to Flicker Elimination Technology and Optix Glare Control, you can count on light that is comfortable for your eyes
- Mini's compact and small design makes it easy to take anywhere, such as a coffee shop or work office
- 90° adjustable stand to provide you with the perfect angle of light

†10,000 lux at 6".

Care and Cleaning

Your HappyLight® Mini is made from high-quality materials designed to last with minimum care. You may want to periodically clean the lamp using a mild, non-abrasive cleaner and soft cloth. When cleaning, make sure you have turned off and unplugged the lamp and allowed sufficient time for it to cool.

Technical Specifications

HappyLight® Mini Cordless LED Light Therapy Light

Dimensions: 5.1" (W) x 5.1" (D) x 0.4" (H)

Environment:

Operating Temperature: -10 to 35°C (14 to 95°F)

Storage Temperature: -20 to 40°C (-4 to 104°F)

Relative Humidity: 10 to 90% non-condensing

Power:

Input Power to the Adapter: 100-240 VAC, 0.5A, 50/60Hz

Input Power to the Lamp: DC5.0V, 3.0A

220V Compatible:

For 220V usage, simply attach a plug adapter.

No converter is needed.

ETL Listed in Canada and the USA

Troubleshooting

If you are experiencing an issue with your HappyLight® Mini, please:

- Make sure the power adapter and plug are inserted fully and securely.
- Make sure there is power to the wall outlet or try another outlet – one that is not connected to a dimmer.

Note: HappyLight® Mini works best with the included cable. Other USB-C cables will not connect completely flush with the USB-C port.

Problem	Check	Solution
Light won't come on.	Have you turned the light on?	The light does not turn on automatically when plugged in. Make sure that you press the power button – (See "How to Use" on page 7).
	Has the plug come loose?	Insert the plug securely into the back of the lamp.
Light Flickers	Has the plug come loose?	Insert plug securely into wall outlet. Insert plug securely into the back of the lamp.
	Is the lamp connected to a dimmer switch?	Connect the lamp plug only to wall outlet without dimmer switch.

Three-Year Limited Warranty

As the leading innovator of light therapy products for over two decades, Verilux, founded in 1956, backs this HappyLight® Energy Lamp with a three-year limited warranty.

Attention: Once opened, please do not return this product to the store where it was purchased for repair or replacement.

Many questions can be answered by visiting www.verilux.com, or you may call our Customer Happiness Team at 1-800-786-6850 during normal business hours.

This limited warranty is provided by Bear Down Brands LLC, 2803 S. Yale St., Santa Ana, CA 92704.

Three-Year Limited Warranty

As the leading innovator of light therapy products for over two decades, Verilux, founded in 1956, backs this HappyLight® Energy Lamp with a three-year limited warranty.

Attention: Once opened, please do not return this product to the store where it was purchased for repair or replacement.

Many questions can be answered by visiting www.verilux.com, or you may call our Customer Happiness Team at 1-800-786-6850 during normal business hours.

This limited warranty is provided by Bear Down Brands LLC, 2803 S. Yale St., Santa Ana, CA 92704.

Verilux warrants this product to be free from defects in material and workmanship for a period of three years from the date of the original retail purchase from Verilux or an authorized Verilux distributor. **Proof of purchase is required for all warranty claims.** During the limited warranty period, Verilux, will, at its option, repair or replace defective parts of this product at no charge to the customer, subject to these limitations: This limited warranty does not include any postage, freight, handling, insurance, or delivery fees. This warranty does not cover damage, defect, or failure caused by or resulting from accident, external destruction, alteration, modification, abuse, misuse, or misapplication of this product.

This warranty does not cover damage to the product resulting from return shipping or handling. Shipping insurance is recommended to help protect your product.

Return Authorization is required for all returns. To obtain a Return Authorization, please contact the Verilux Customer Happiness Team.

If, during the three years of ownership, this product fails to operate properly, it should be returned as specified at verilux.com/warrantyreplacement or as instructed by a Verilux Customer Happiness Specialist.

Due to ongoing improvements, actual product may have slight variations from the product described in this manual.

Please visit our website at: verilux.com
or contact our Customer Happiness Team at

1-800-786-6850 hello@verilux.com

Specialists are available Monday - Friday, 8am-5pm PST

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La garantie limitée est offerte par Bear Down Brands LLC, 2803 S. Yale St., Santa Ana, CA 92704.